



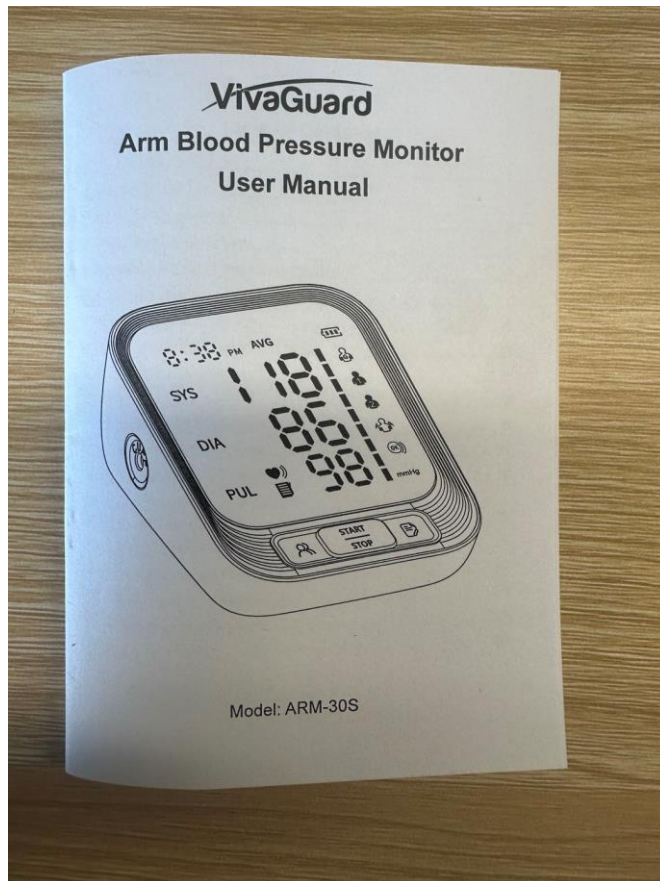


Table of Contents	
1. Unpacking Inspection.....	1
2. Packing List.....	1
3. Safety Precautions	2
4. Product Composition.....	3
5. Intended Use / Instructions for Use	3
6. Contraindication	3
7. Product Parts.....	4
8. WHO Blood Pressure Indicator.....	5
9. Power Connection	5
10. Function Setting.....	6
11. How to Apply The Arm Cuff	9
12. How to Take Proper Measurements	11
13. Specifications	15
14. Contraindications, Precautions, Warnings and Prompt Instructions.....	17
15. Common Q & A on Blood Pressure.....	20
16. Abnormal Phenomena and Handling	22
17. Cleaning and Disinfection	24
18. Upkeep and Maintenance.....	26
19. Appendix 1 EMC Information.....	27

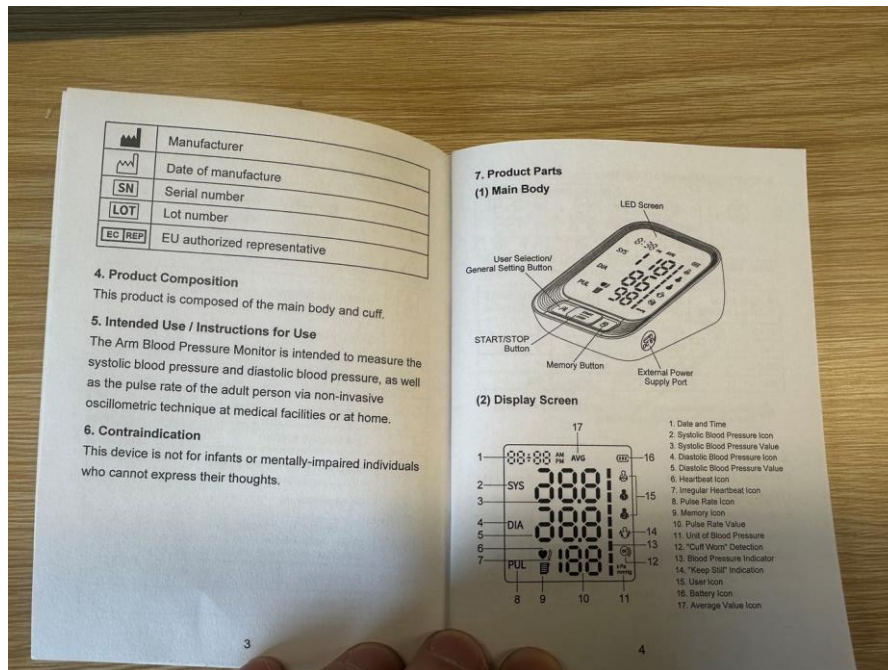

8. WHO Blood Pressure Indicator

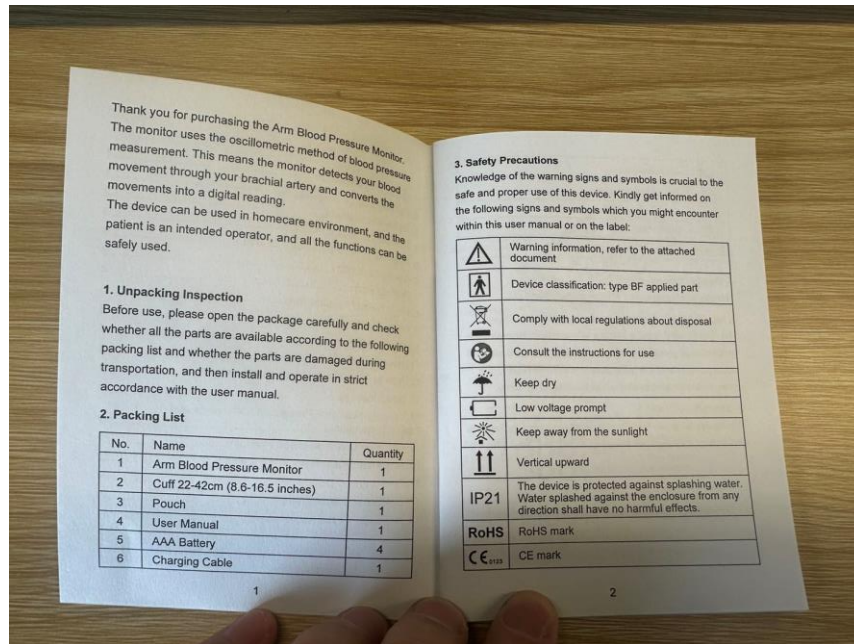

Systolic Blood Pressure (mmHg)	Diastolic Blood Pressure (mmHg)	Color Display	Hierarchical Relationship
≥160	≥100	red	or
140-159	90-99	yellow	or
90-139	60-89	green	or
<90	<60	yellow	and

Warning: Never diagnose or treat yourself based on the readings. Please always consult with your physician.

9. Power Connection

(1) Install The Batteries

- 1) Open the battery cover according to the method shown in the figure.
- 2) Place 4 AAA batteries in the battery compartment, and pay attention to the electrode indication of the batteries. Install the battery as indicated in the picture.



5

Remove the batteries if you do not intend to use the device for a long time (over 3 months).

(2) Type-C Connection for Power Supply


In addition to the batteries, the power can also be supplied by plugging into d.c. 5V external power supply through a Type-C port. (The power supply line is not included in the packing list.)



10. Function Setting

(1) User Selection


In the power-off mode, press the "M" button to enter the user group selection interface. Then press the "M" button again to switch and select user groups.

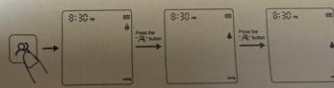

6

(2) Date And Time Setting

In the power-off mode, press the "M" button for about 3 seconds to enter the date setting interface, and the "year" will flash. Press the "M" button to adjust the year, and press the "M" button to confirm the selection. When the year is set, it will automatically enter the month setting. Press the "M" button to adjust the month, and press the "M" button to confirm the selection. Follow the same steps to adjust the date/hour/minute.

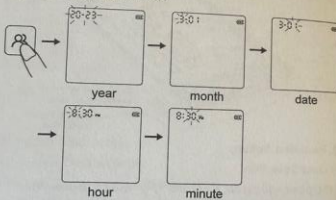


7



(3) Unit Display Setting

There are two units of blood pressure display, mmHg and kPa. The default unit is mmHg.

In the power-off mode, press the "START/STOP" button for about 5 seconds to enter the unit selection. Press the "M" button to switch between mmHg and kPa, and then press the "M" button to confirm the selection.

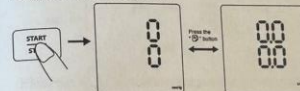


8

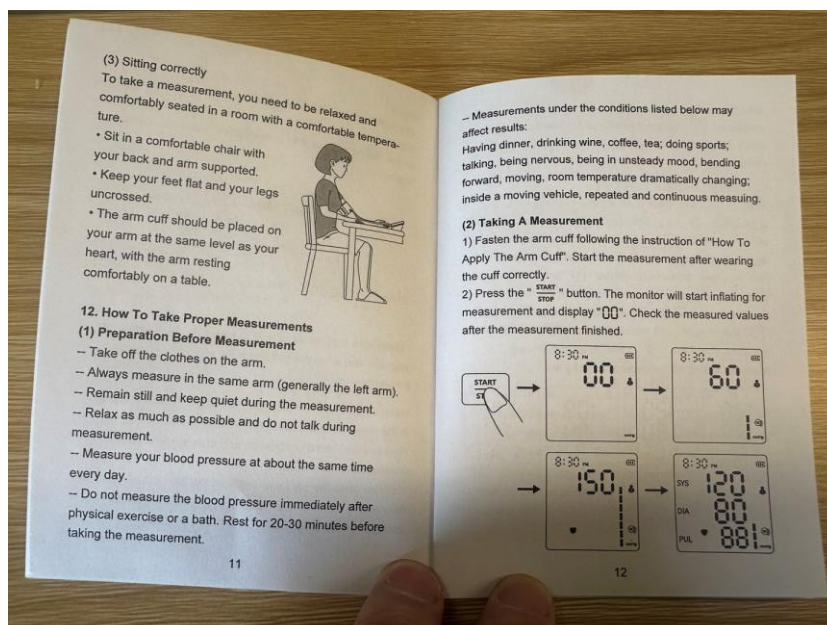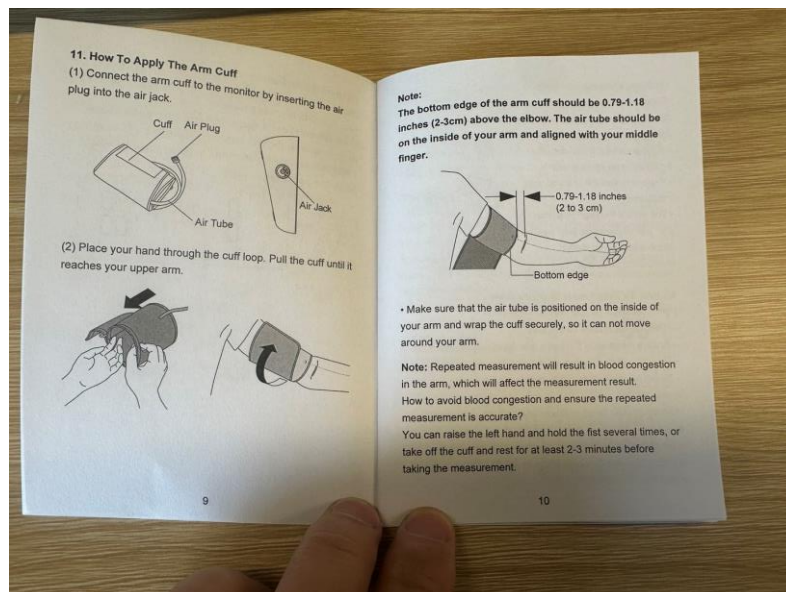

Note: If you feel uncomfortable during the measurement, press the "START/STOP" button immediately to stop the measurement. Please consult your doctor if unexpected readings are obtained.

(3) Memory Function

- 1) Each measured value is stored automatically under the appropriate user group. This device can store up to 199 sets of measurements for user 1 and user 2. Once the memory log is full, old values will be refreshed with new ones.
- 2) In the power-off mode, press the "M" button once and the device will display the average value of the blood pressure measurements of the last 2 or 3 times. Press the "M" button again, and the latest measured value will be displayed. Press the "M" button again and the rest measured values will be displayed one by one.

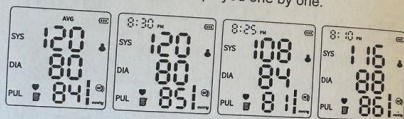

13

3) Delete Memory

In the power-off mode, press the "M" button to select the user group whose measured values need to be deleted. Press the "START/STOP" button to power off the device and press the "M" button once to activate the screen. Then keep pressing the "M" button for about 3 seconds to delete the memories of the selected user and the "no" icon will appear on the screen.



(4) "Cuff Worn" Detection

The "OK" icon is always displayed on the screen when the cuff is wrapped correctly. When the cuff is too loose, the "OK" icon will always flash to remind you. If the "OK" icon is flashing, please press the "START/STOP" button to stop the measurement.

14

(5) "Keep Still" Indication

The "Keep Still" icon flashes when you move your body or shake your arm during the measurement, which may cause incorrect measurement results. Please adjust your posture and measure again.



13. Specifications



Model	ARM-30S	
Display	LED screen	
Measuring Method	Oscillometric measurement	
Measuring Part	Upper arm	
Pneumatic Pressure Measuring Range	0-295 mmHg (0-39.3 kPa)	
Maximum Pressure Protection	295 mmHg (39.3 kPa)	
Measurement Range	Blood pressure value	SYS: 57-255 mmHg (7.6-33.4 kPa); DIA: 25-195 mmHg (3.33-26 kPa);
	Pulse rate	40-199 bpm
Accuracy	Blood pressure value	±3 mmHg (±0.4 kPa)
	Pulse rate	±5%
Low Battery	When the power is lower than 4V±0.1V, the device will be turned off automatically.	
Power Source	4*AAA batteries or d.c. 5V power supply line	

15

Memory	2 users x 199 memories + guest mode	
Dimension	130 mm (L) x 100 mm (W) x 49 mm (H) (5.12 inches x 3.94 inches x 1.93 inches)	
Screen Size	67 mm (L) x 68 mm (W) (2.6 inches x 2.7 inches)	
Cuff Size	22-42 cm (8.66-16.53 inches)	
Weight	About 221g (without batteries)	
Anti Electronic Shock Type	Internal Power Supply	
Auto Power-off	1 minute without operation	
Anti Electronic Shock Degree	Type BF	
Operation Mode	Continuous operation	
Protection Against Harmful Ingress of Water or Particular Matter	IP21	
Service life	5 years	
Protection Against Electric Shock	Internally powered supply	
Operating Environment	Temperature condition	5°C-40°C
	Humidity condition	15%-90%RH
	Atmospheric condition	70kPa-106kPa

16

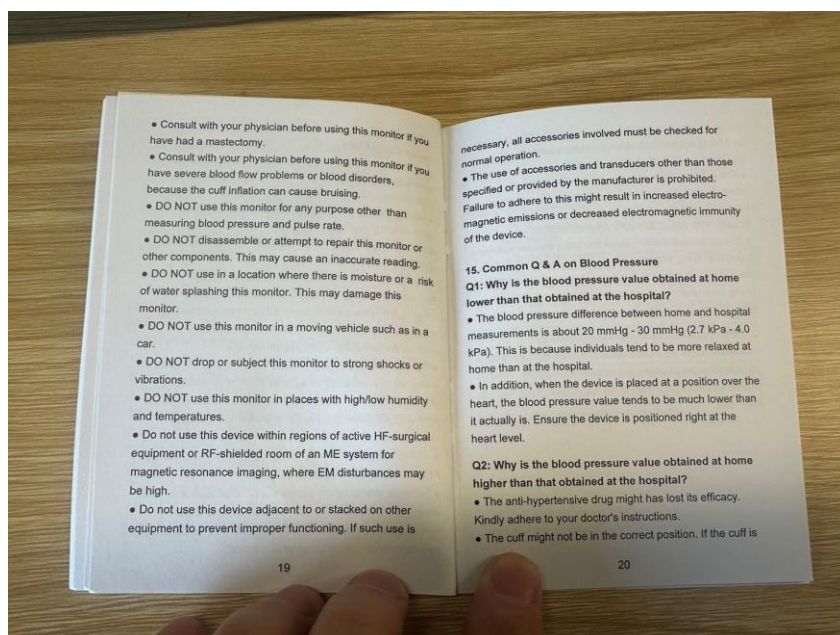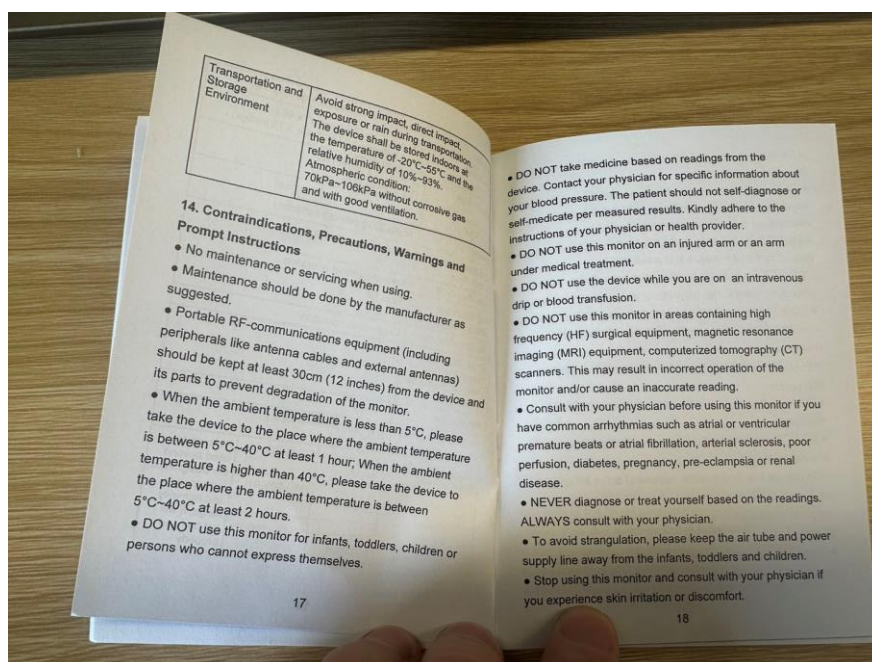

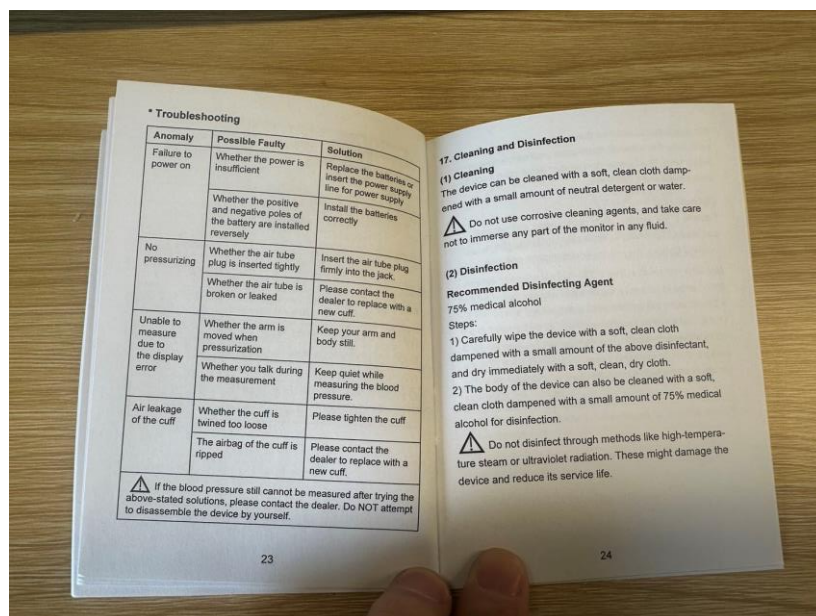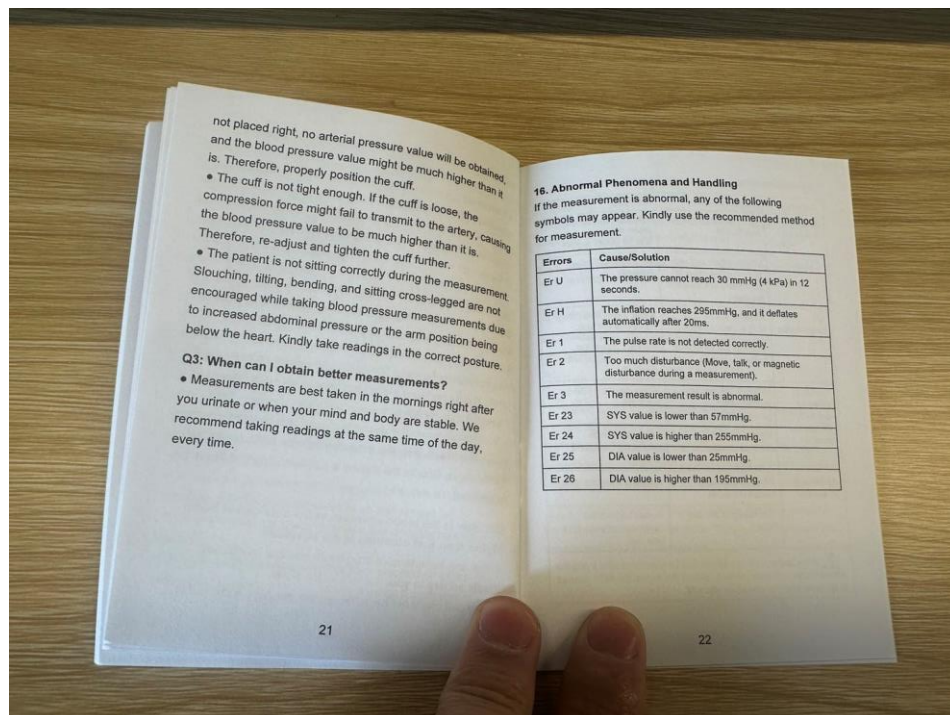

It is suggested to disinfect the monitor before and after use each time. Each time of disinfection shall be completed within 1min. The number of repeated disinfection each time shall not exceed 2 times.

(3) Disposal

Dispose of the monitor, other components and optional accessories according to applicable local regulations. Unlawful disposal may cause environmental pollution.

Notes

- Do not bend or crease the air tube excessively.
- Do not store the monitor or its components:
 - if the monitor or its parts is wet.
 - in locations with extreme temperatures, humidity, direct sunlight, dust, or corrosive gases.
 - in areas with a high risk of vibrations or shocks.

18. Upkeep and Maintenance



Water or neutral detergent



- Always keep the surface of the device clean and tidy, helpful to prolong its service life.
- If the device is dirty, please wipe with a dry soft cloth. If the dirt cannot be eliminated easily, wipe with a soft cloth stained with water or neutral detergent, and then dry with a dry cloth.
- No maintenance or maintenance required when operating the device.

Do not allow water or other liquids to flow into the device.

19. Appendix 1 EMC Information

Guidance and manufacturer's declaration - Electromagnetic emission



The Arm Blood Pressure Monitor is intended for use in the electromagnetic environment specified below. The customer or the user of the Arm Blood Pressure Monitor should assure that it is used in such an environment.

Emissions	Compliance	Electromagnetic environment - guidance
RF emissions CISPR 11	Group 1	The Arm Blood Pressure Monitor uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class B	The Arm Blood Pressure Monitor is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
Harmonic emissions IEC61000-3-2	N.A.	
Voltage fluctuations/flicker emissions IEC61000-3-3	N.A.	

Guidance and manufacturer's declaration - Electromagnetic immunity

The Arm Blood Pressure Monitor is intended for use in the electromagnetic environment specified below. The customer or the user of the Arm Blood Pressure Monitor should assure that it is used in such an environment.

Immunity test	IEC 60601 test level	Compliance level
Electrostatic discharge (ESD) IEC 61000-4-2	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air
Electrostatic discharge (ESD) IEC 61000-4-2	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air
Electrical fast transient/burst IEC 61000-4-4	Not applicable	Not applicable
Surge IEC 61000-4-5	Not applicable	Not applicable
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	Not applicable	Not applicable
Power frequency Magnetic field IEC 61000-4-8	30 A/m, 50/60Hz	30 A/m, 50/60Hz
Conducted RF IEC61000-4-6	Not applicable	Not applicable

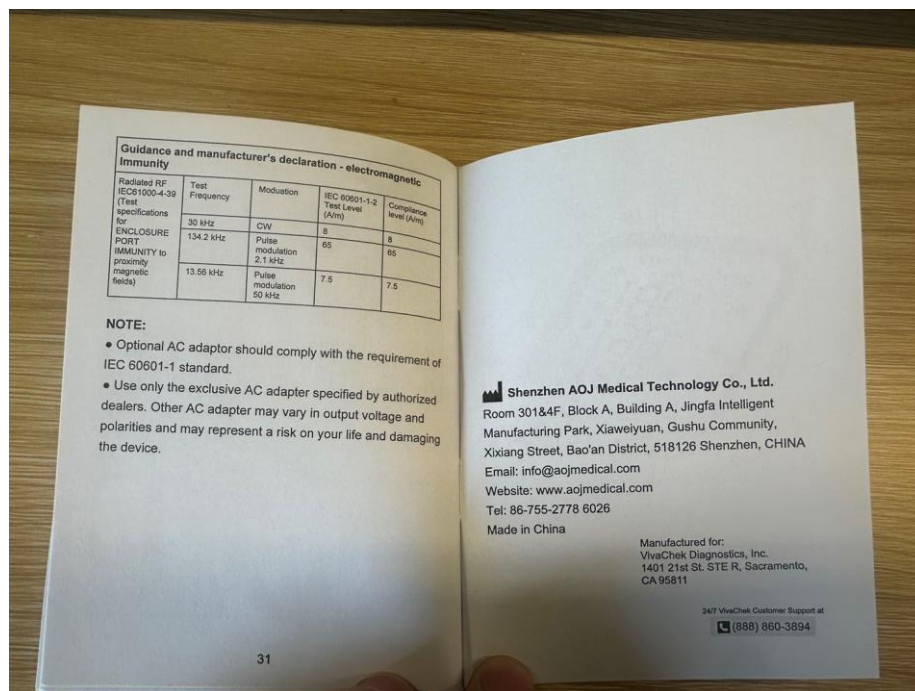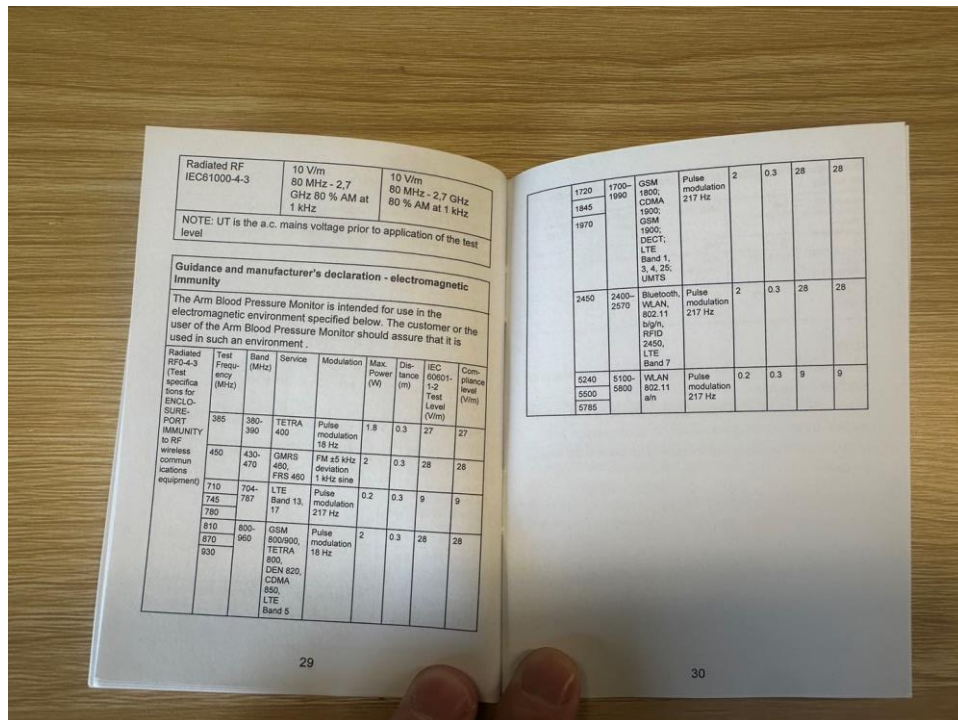



PN: 1201046701
Manual version: A0
Software version: V1.0
Revision date: 2022-08